

# HU-FRIEDY BAGETTE

Self-Sealing Sterilization Pouches

Hu-Friedy's Bagette self-sealing sterilization pouch line includes a variety of sizes created to fit all sterilization needs:

- Four pouch sizes specifically designed to fit all of Hu-Friedy's IMS Signature Series cassettes.
- Additional sizes are available for the sterilization of other items including rubber dam clamps, burs, handpieces and loose instruments.



## POINTS OF PERFORMANCE

Internal and external indicators for both Steam and EO Sterilization. Responds to all variables of sterilization: time, temperature & presence of steam

Thumb notch for ease of opening

High quality medical grade paper and multi-layer blue tinted film provide strength and durability

Thick, strong thermal seal to protect against cassette and instrument breaches

Wide self-seal adhesive strip to provide reliable seal

Pouches are lead free and are not made with natural rubber latex

**FOUR NEW BAGETTE SIZES  
TO CHOOSE FROM!**

How the best perform




# COLOR CHANGE GUIDE


Below is a color change guide to aid in the interpretation of the color change that occurs when the pouches are processed in an autoclave, or ethylene oxide sterilization chamber.

Ideally, the indicator should change color from the first box to the last box, but due to the large variation in sterilization chambers, this is not always the case. Therefore, there is a suggested acceptable range of color change. See diagrams below for details.

### STEAM STERILIZATION

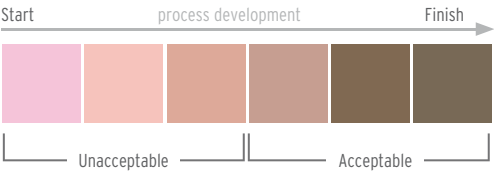

### EO STERILIZATION

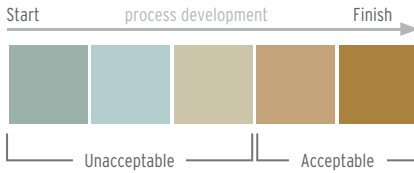

Caution: the chemical indicators are to be used as a sterilization process indicator only. The indicators should be used in conjunction with biological indicators (spore strips) to ensure the devices within the pouches are sterile.

## BAGETTE USAGE CHART

PART CODE	SIZE	POUCHES/BOX	SUITABLE CONTENTS
IMS-1346	2-1/4" x 4" (57 mm x 102 mm)	200	Rubber dam clamps or burs
IMS-1347	2-3/4" x 9" (70 mm x 229 mm)	200	1 handpiece or 1-3 hand instruments
IMS-1236	3-1/2" x 9" (89 mm x 229 mm)	200	2-5 hand instruments
IMS-1348	5-1/4" x 10" (133 mm x 254 mm)	200	Handpiece or multiple instruments
IMS-1237	7-1/2" x 13" (190 mm x 330 mm)	100	<ul style="list-style-type: none"> <li>• 1/4 DIN Cassette – 5 instruments; IME14DIN5X</li> <li>• Exam Cassette – 5-instrument; IM605X</li> <li>• Double-Decker Cassette – 8 instrument; IM908X</li> <li>• Large Handpiece Cassette; IM60HX</li> <li>• Small Ultrasonic Insert; IM30HUX</li> <li>• Symmetry IQ® Handpiece Cassette; IM60IQ8, IM30IQ8</li> </ul>
IMS-1238	10" x 14" (254 mm x 356 mm)	100	<ul style="list-style-type: none"> <li>• 1/2 DIN Cassette – 8-10 instruments; IME12DINXXX</li> <li>• Small Cassette – 8 instruments; IM508X</li> <li>• Small Cassette – 10 instruments; IM510X</li> <li>• Double-Decker Cassette – 14 instrument; IM914X</li> <li>• Radiography Cassette; IM7000</li> <li>• Space Saver Seven Cassette; IMSS7X</li> <li>• Small Thin Ortho Cassette; IMORTHSMB</li> <li>• Small Resin Cassette – 7-16 instruments; IMS-3107, IMS-3112, IMS-3116</li> </ul>
IMS-1239	12" x 17" (305 mm x 432 mm)	100	<ul style="list-style-type: none"> <li>• DIN Cassette – 13-20 instruments; IMEDINXXX</li> <li>• Large Cassette – 16 instruments; IM416X</li> <li>• Large Cassette – 20 instruments; IM420X</li> <li>• Double-Decker Cassette – 20-instrument; IM920X</li> <li>• Signa-Stat Cassette – 12-instrument; IM812X</li> <li>• Large Thin Ortho Cassette; IMORTHLG8</li> <li>• Large Resin Cassettes – 12-24 instruments; IMS-1118, IMS-1124, IMS-1112, IMS-1116</li> </ul>
IMS-1345	13" x 20" (330 mm x 508 mm)	100	<ul style="list-style-type: none"> <li>• Oral Surgery Cassette; IMOSX</li> <li>• Large Custom Cassette; I18</li> </ul>

Visit us online at [Hu-Friedy.eu](http://Hu-Friedy.eu)

Customer Care Department: Free Phone 00800 48 37 43 39  
Free Fax 00800 48 37 43 40 | Email: [info@hu-friedy.eu](mailto:info@hu-friedy.eu)

©2014 Hu-Friedy Mfg. Co., LLC. All rights reserved. IMS-520HFE/042014

